

PORTABLE MESH NEBULIZER

MODEL NO: HNK-MESH-01

MADE IN
TURKEY



nimo®



NIMO® Portable Mesh Nebulizer is designed for the treatment of asthma, allergies and other respiratory ailments. Thanks to the active mesh technology that converts electrical energy into mechanical energy, HNK-MESH-01 transforms liquid drugs into aerosol mist and allows them to be taken into the body by respiration. NIMO® Portable Mesh Nebulizer, which has titanium alloy active mesh technology, can give the drug to the patient quietly and quickly without disturbing the patient. The device is suitable for use by patients of all ages. HNK-MESH-01 is easy to use. HNK-Mesh-01 is lightweight and ergonomically designed for easy use at home, at work or on travel.

A companion which is with you whenever you want

High
Density
Steam
Output

Low
Power
Consumption
and Quiet
Operation

60 min
Usage With
Full Charge

Rechargeable
Li-Titanate
Battery

Full Charge
in 90 Minutes
With USB
Cable

Lightweight
and Portable
Mini Compact
Design

Automatic
Shutdown in
10 Minutes

Easily
Attachable,
Cleanable
Medicine
Container

Minimum
Drug Residue



Fast Treatment Opportunity Thanks to High Nebulization Rate

Thanks to the mesh technology with 2480 grain pores in 2.5 micron size, the vaporized drug can be inhaled by the parts of the lungs where it will be most effective

USB Cable Charging Feature

Thanks to its rechargeable battery, it puts an end to the need for sockets. After charging the battery of this low-power and extremely quiet device, you can have it with you wherever you want.





Easy-to-Use Medicine Container

In addition to easy attachable design, it will provide you comfort with its easy cleaning

We Care Your Time

It offers you ease of use with its lightweight and portable design in travel situations. You can take the NIMO® Portable Mesh Nebulizer with you wherever you go



Medicine
Container Lid

Titanium
Alloy Mesh

Medicine
Container

Medicine
Container
Eject Button

The Main Body
of Device

On/Off
Button

USB Charge
Port

Adult and
Paediatric
Mask



Nozzle



Nimo® Portable Mesh Nebulizer, Your Companion in Fast and Efficient Treatment

Thanks to its light and ergonomic structure, it does not disturb the patient. You can get comfortable by using this device, which has many advantages such as being able to give the medicine to the patient silently and quickly.

Led Indication Light

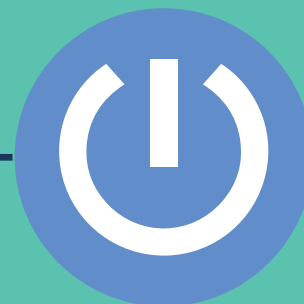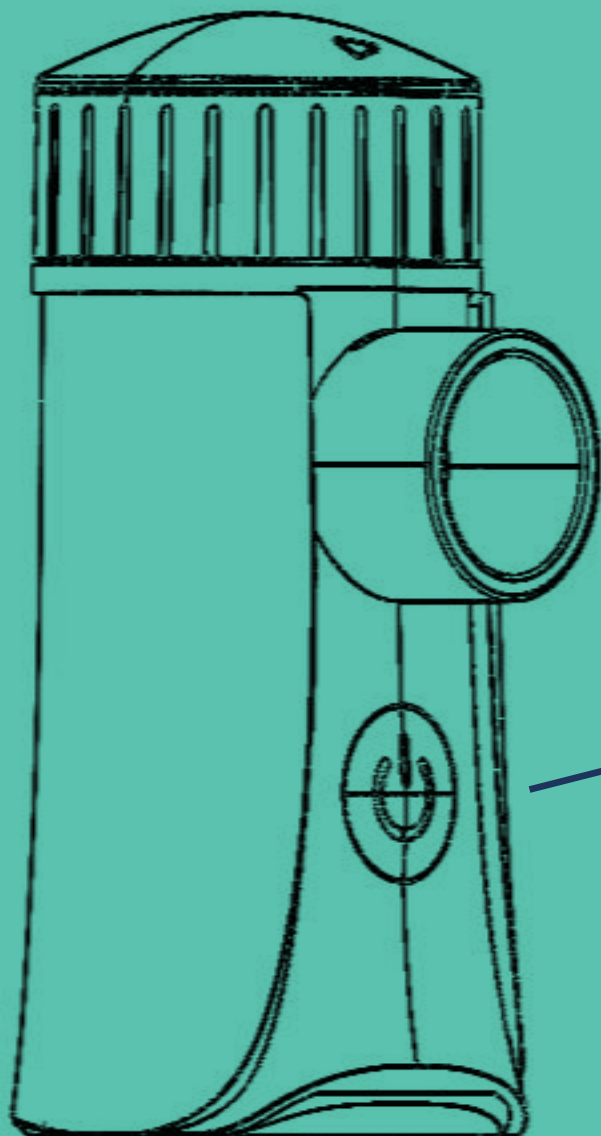

Blue Light:
Device is Working



Green Light:
Device is Charging



No lights:
Device Closed

TECHNIC FEATURES

Working Method	Active Vibrating Mesh Technology
Rated Operating Voltage	DC 4.8V (built-in lithium titanate battery) or DC 5.0V 1A
Power Consumption	<4.0 W
Nebulisation Speed	0.5 ± 0.15 ml/min
MMAD	2.5±30%µm
Vibration Frequency	120 kHz±10 kHz
Medicine Container Capacity	10 ml
Type	Type BF Applicable Parts
Solid and Liquid Protection Class	IP22
Size	41mm(L)x45mm(W)x91mm(H)
Weight	110 gr (Including Battery)



BOX CONTENTS

1 NIMO Mesh Nebulizer	1 Piece Nozzle	1 Piece Adult Mask	1 Piece Child Mask
1 Piece Micro USB Cable	User Manual	Warranty Card	1 Piece Transport Bag



T.C. SAĞLIK BAKANLIĞI
Onaylı Ürün



ÜRÜN TAKİP SİSTEMİ



All the best
FOR YOUR HEALTH!